Oral-B®

ProfessionalCare™

8000 series






powered
by **BRAUN**

Type 4729
Modelo 4729
Modèle 4729

3-728-151/03/II-06
USA/CDN

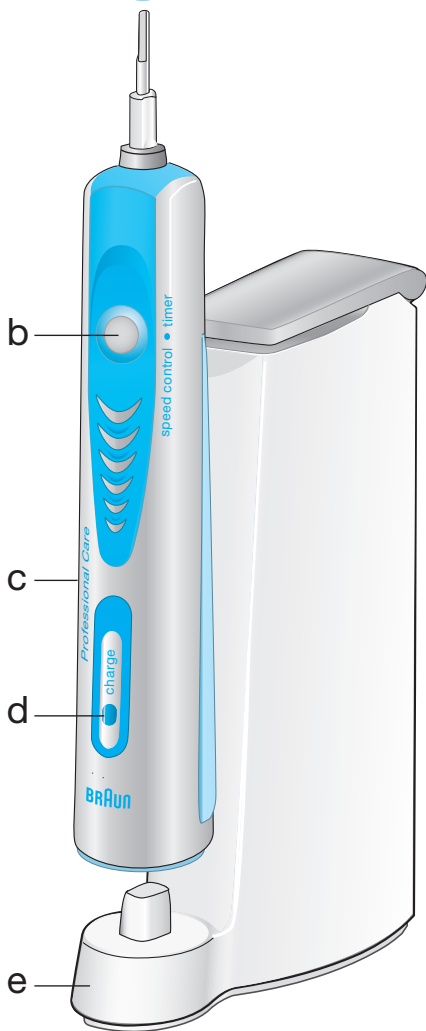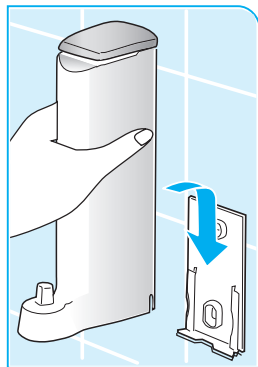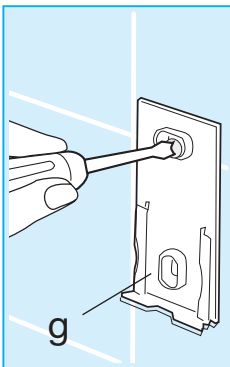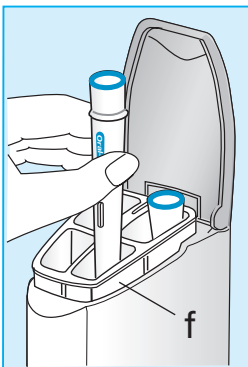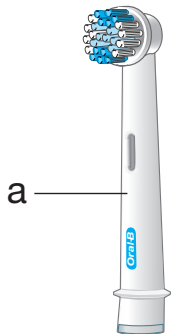
h

6

6

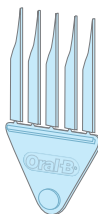6

Table of contents / Tabla de contenidos / Table des matières

English	5
Important Safeguards	6
Introducing your Oral-B® ProfessionalCare™ 8000 series	7
Connecting and charging	7
Using your toothbrush	8
Individual brushheads and accessories	9
Cleaning recommendations	11
Environmental notice	11
Questions? Ask Oral-B	12
Limited 2-Year Warranty	13
The 30-Day Challenge	15
Español	16
Medidas de Seguridad	17
Conociendo a su Oral-B® ProfessionalCare™ 8000 series	18
Conexión y carga	18
Usando su cepillo	19
Cabezales individuales y accesorios	20
Recomendaciones de limpieza	22
Nota respecto al medio ambiente	22
¿Preguntas? Hagaselas a Oral-B	23
Garantía de 2 años	23
El reto de 30 días	24
Français	25
Importantes mesures de sécurité	26
Renseignements sur la Oral-B® ProfessionalCare^{MC} série 8000	27
Branchement et charge	27
Utilisation de la brosse à dents	28
Têtes de brosse et accessoires particuliers	29
Recommandations de nettoyage	31
Avis environnemental	31
Des questions ? Demandez à Oral-B	32
Garantie restreinte de deux ans	33
Le défi de 30 jours	34

Congratulations!

You have chosen the best power toothbrush technology available.

Oral-B® ProfessionalCare™ 8000 series has been developed together with leading dental professionals as a complete oral care system for superior cleaning, natural whitening, gum care and freshening of your whole mouth.

It reverses gingivitis for firm and healthy gums, cleans below the gumline and behind back teeth, polishes away stains and freshens breath for a healthier, brighter smile.

It delivers the ultimate in oral care through superior performance:

- Patented 3D oscillating-pulsating technology – independently acclaimed to be the superior way to clean teeth and care for gums
- Innovative accessories and brushheads to clean and polish teeth as well as freshen your tongue and mouth
- Customized Speed Control which offers a broad variety of speed levels to meet your individual brushing needs
- Premium handle design and ergonomics for greater comfort and control

Oral-B ProfessionalCare 8000 series is brought to you by Oral-B, the worldwide leader in toothbrushes and the brand more dentists use themselves.

If you have any questions about using your new **Oral-B ProfessionalCare 8000 series**, please visit our web site at www.oralb.com or call our toll free telephone number: 1-800-566-7252

Thank you for choosing
Oral-B ProfessionalCare 8000 series.

Sincerely,

Bruce Cleverly
President, Oral-B

Oral-B® ProfessionalCare™ 8000 series.
Heals, restores, renews.™

IMPORTANT SAFEGUARDS

Oral-B ProfessionalCare 8000 series has been carefully designed to offer you and your family a new and unique brushing experience that is both safe and effective.

However, when using electrical products, basic safety precautions should always be followed. **Please read carefully before use.**

DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop or insert any object into any opening.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charger has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.

SAVE THESE INSTRUCTIONS

Introducing your Oral-B® ProfessionalCare™ 8000 series

- a FlexiSoft® brushhead
- b On/off button
- c Oral-B ProfessionalCare handle
- d Charge indicator
- e Charging unit
- f Brushhead compartment
- g Wall holder
- h Range of brushheads and accessories


Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and can be used in the bathroom without hesitation.


- Plug the charging unit (e) into an electrical outlet and place the toothbrush handle (c) on the charging unit.
- The blue charge indicator (d) flashes while the handle is being charged. If you then use your brush and place it back on the charging unit, the charging light will again flash to let you know it is charging.
- When the handle is fully charged, the charge indicator shines continuously. A full charge takes at least 16 hours and enables approximately two weeks of regular brushing.
- When the rechargeable battery is running low, your toothbrush will slow down, then stop.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle by regular use at least every 6 months.

Specifications

Refer to the bottom of the charging unit for voltage specifications.

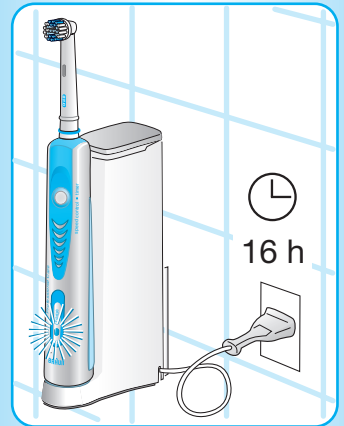Using your toothbrush


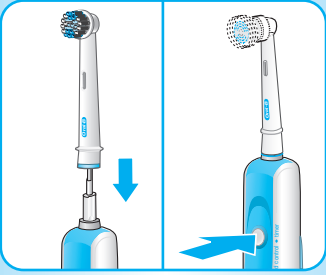
Attach one of the brushheads (1–3) to the handle. They can be used with any kind of toothpaste. To avoid splashing, guide the brushhead to your teeth before switching the appliance on.

Customized Speed Control


Oral-B ProfessionalCare 8000 series offers a variable range of speeds to best meet your individual needs. High speed is ideal for tooth and gumline cleaning. For cleaning sensitive teeth and gum tissue as well as precision use of the accessories (4 and 5), we recommend using a lower speed.

- Press the on/off button once to operate at maximum speed.
- Continue to hold the on/off button to shift from the highest speed to the Customized Speed Control mode.
- As long as you keep the on/off button pressed, the handle will continuously run through all speed levels.
- Release the on/off button to brush at your preferred speed level.
- Push the on/off button once more to turn the toothbrush off.
- If you then turn it on again you will automatically come back to the maximum speed level.

outsides

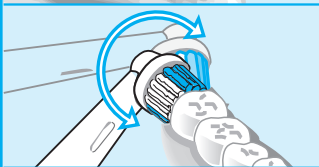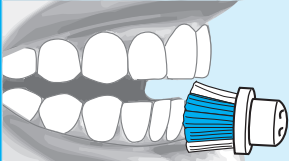insides

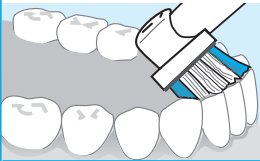tops

Ultimate Plaque Removal

For cleaner teeth and healthier, firmer gums, **Oral-B ProfessionalCare 8000 series** gently cleans in areas where plaque bacteria can collect. It removes up to 97 % of plaque from hard-to-reach areas, prevents build-up of tartar and is gentle on teeth and gums.

Brushing technique:

- Guide the brushhead slowly from tooth to tooth.
- Hold the brushhead in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

Pressure Sensor

For optimal brushing results **Oral-B ProfessionalCare 8000 series** has a pressure control feature to ensure that the right degree of pressure is applied. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brushhead will continue but the pulsation will stop. In addition to feeling it, you can also hear a different sound.

Professional Timer

A short stuttering sound at 30-second intervals, which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a better overall clean.

2-Minute Timer

A long stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.

During the first few days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist.

Individual brushheads and accessories

The **Oral-B ProfessionalCare 8000 series** offers you a variety of brushheads and accessories, allowing you to select those which best meet your individual needs. They are all designed for daily use:




1 FlexiSoft® brushhead


The FlexiSoft® brushhead features the superior oscillating and pulsating action to effectively clean your teeth, clean along the gumline and in hard-to-reach areas such as behind teeth. The compact round brushhead surrounds each tooth for a complete clean. The green FlexiSoft® bristles bend, enabling both gentle brushing and deeper reach of the blue Interdental Tips between teeth. Apply the brushing technique as described in the section «Ultimate Plaque Removal».

2 Oral-B PowerPolisher™ brushhead (not with all models)

The PowerPolisher brushhead is specifically designed to polish and whiten your teeth by removing stains from coffee, tea and tobacco. It features the superior oscillating and pulsating head action of the FlexiSoft® brushhead to effectively clean teeth, combined with the special polishing cup for effective brightening.

Apply the same brushing technique as described in the section «Ultimate Plaque Removal». We do not recommend using the PowerPolisher brushhead with braces.

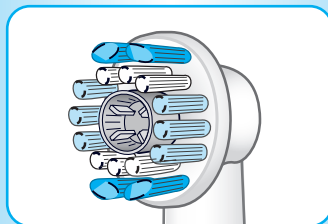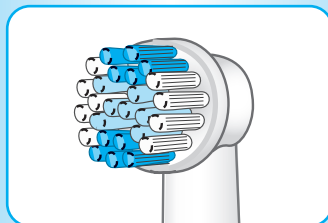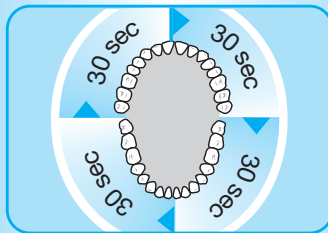
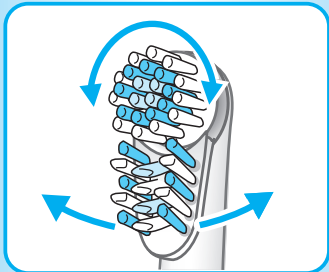
3 Oral-B DualAction™ brushhead (not with all models)

The DualAction brushhead offers two types of cleaning motion in one. It combines the superior oscillating and pulsating motion with an additional side-to-side motion, for more cleaning action. The oscillating motion stimulates gums while removing plaque and bacteria. Multi-angled bristles clean between teeth and below the gumline. The larger brushhead results in increased freshening sensation.

With this brushhead you can either move from tooth to tooth or apply the same brushing technique as with any manual toothbrush.

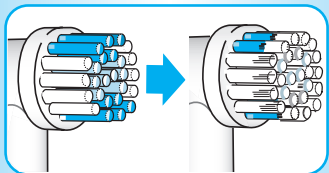
All three brushheads feature blue INDICATOR® bristles to help you monitor brushhead replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway within 3 months, signaling it is time to change the brushhead.

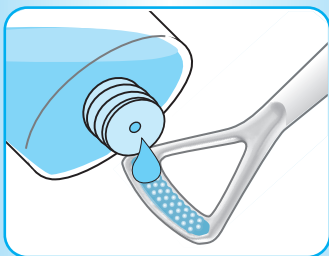4 Tongue Freshener (not with all models)

The Tongue Freshener freshens breath by removing odor-causing bacteria from the tongue. It can be used everyday – before or after brushing. You can use mouthwash with the Tongue Freshener.


Attach it to the handle, turn the appliance on and drag it across the tongue from back to front. We recommend a low speed setting for enhanced control during use. Rinse after use. Replace the Tongue Freshener when the material becomes worn.

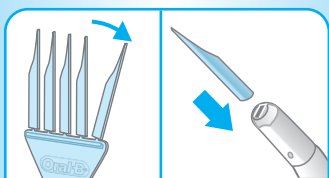
5 Interdental Cleaner (not with all models)

The Interdental Cleaner effectively cleans between the teeth, protects against gum disease and reduces gingivitis. It should be used before brushing.

Break one green cleaner off the base (6) and insert it into the slot in the neck of the white pick holder. Attach it to the handle and turn the appliance on. Depending on which part of your mouth you are cleaning, the pick holder should be adjusted by turning it, so that the blunt side always faces your gums. Guide the cleaner between your teeth and clean systematically from space to space.

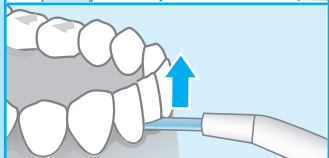
We recommend a low speed setting for enhanced control during use. The Interdental Cleaner can be used approximately two to three times until it becomes too soft to effectively fit and clean between the teeth.

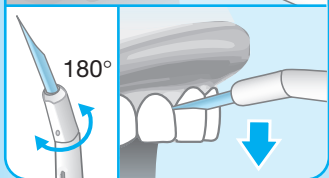Oral-B ProfessionalCare 8000 series replacement brushheads and accessories are available in stores that sell Oral-B products, at the Essentials.com web site and at Braun Service Centers.

Cleaning recommendations

- After brushing, rinse your brushhead or accessory thoroughly under running water for several seconds with the handle switched on.
- Switch off the handle and remove the brushhead or accessory from the handle. Clean all parts separately under running water, then wipe them dry before you reassemble the brush and put the handle on the charging unit.
- From time to time, also clean the charging unit and the brushhead compartment with a damp cloth.
- For cleaning purposes, the brushhead compartment is removable.

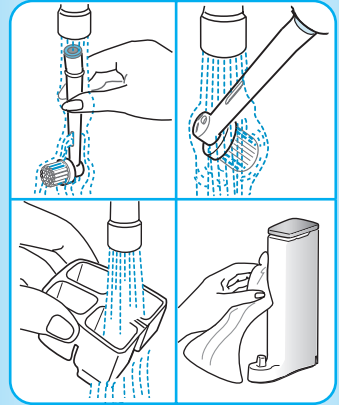Environmental notice


This product is provided with a nickel-hydrate battery, which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Center or at approved recycling or disposing locations in accordance with local regulations.

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution! Opening the handle will destroy the appliance and invalidate the warranty.

Subject to change without notice.

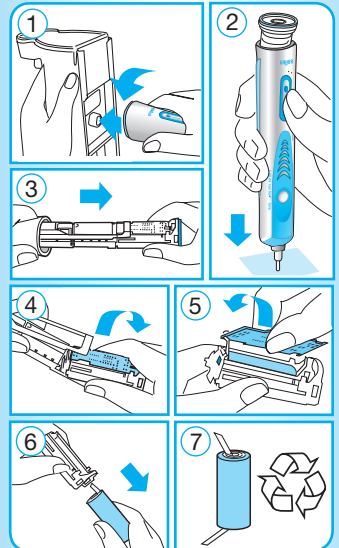Questions? Ask Oral-B.

I have sensitive teeth and gums. What is the best way to brush with my Oral-B ProfessionalCare 8000 series?

Oral-B ProfessionalCare 8000 series is as safe as a soft manual toothbrush at all speed levels. However, it gives you the flexibility to choose your desired speed level, to best meet your individual needs.

My dentist told me I have gingivitis. Can my Oral-B ProfessionalCare 8000 series help?

Yes. With its patented 3D cleaning action **Oral-B ProfessionalCare 8000 series** has been clinically proven to prevent and even reverse gingivitis, an early form of gum disease, and to remove twice as much plaque as a manual toothbrush. It thereby delivers superior results versus manual brushing with minimum effort.

I wear braces. What is the recommended way to brush?

Gently place the brush at an angle between the brackets and gums. Using a slower speed level will help maximize the precision. The Tongue Freshener and Interdental Cleaner can be used with braces, but use of the PowerPolisher brushhead is not recommended. In addition, Oral-B's Ortho brushhead is specifically designed to clean around braces and wires and can be used with the **Oral-B ProfessionalCare 8000 series**.

I have bridges, crowns and fillings in my mouth. What is the best way to brush with my Oral-B ProfessionalCare 8000 series?

Like all Oral-B power toothbrushes, **Oral-B ProfessionalCare 8000 series** can be used by consumers with bridges, implants and crowns. We recommend using a slower mode for even better control.

I have recently had oral surgery. Can I use Oral-B ProfessionalCare 8000 series?

You should always follow the instructions given by your oral surgeon or periodontist after an operation. Typically in the first few days following surgery, teeth and gums may be too sensitive for any toothbrush use. As quickly as your teeth and gums heal to permit brushing, you can proceed with using **Oral-B ProfessionalCare 8000 series**. Its slow mode of brushing speed may be beneficial. The superior plaque removal capability will be particularly valuable, as a clean mouth will aid the healing process.

What happens if I use the Oral-B ProfessionalCare 8000 series Interdental Cleaner and Tongue Freshener at a high speed instead of the recommended low speed?

The **Oral-B ProfessionalCare 8000 series** accessories can be operated at high or low speeds. However, we recommend the lower speed for optimal control.

If I have comments/questions about my Oral-B ProfessionalCare 8000 series, what should I do?

For any additional questions or comments you have, please visit our website at www.oralb.com or call our toll free number: 1-800-566-7252

FOR USA ONLY

Limited 2-Year Warranty

Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts and accessories:

Contact an authorized Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To order original Braun replacement parts and accessories (including replacement brushheads) by phone:

Please call 1-800-566-7252.

To speak to an Oral-B Consumer Service Representative:

Please call 1-800-566-7252.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service centre, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts & accessories:

Contact an authorized Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To speak to an Oral-B Service Representative:
Please call 1-800-566-7252.

The 30-Day Challenge

Take the challenge! Experience the **Oral-B ProfessionalCare 8000 series** difference!

Oral-B ProfessionalCare 8000 series is more than a toothbrush. It is a complete oral care system for superior cleaning, natural whitening, gum care and freshening of your whole mouth.

No wonder that Oral-B is the brand more dentists use themselves worldwide. In fact, we are so convinced you'll fall in love with your new **Oral-B ProfessionalCare 8000 series** that we give you this risk-free offer.

Try **Oral-B ProfessionalCare 8000 series** for 30 days. If you are not 100% satisfied, receive a full refund. No questions asked. Guaranteed.

For a full refund, please return your handle, refills, and charging unit with your original cash register receipt. For further details, please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

U.S. only:

Attention: 30 Day Challenge
Braun Division
The Gillette Company
Sobin Park and A Street, MS 1F
South Boston, MA 02127

Canada only:

Attention: 30 Day Challenge
Braun Products
Gillette Canada Company
4 Robert Speck Parkway
Mississauga ON L4Z 4C5

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